



Yes, you can.[®]

Invacare[®] Action³

Roll into action

The *Invacare[®] Action³* has been designed to enable easier assessment, and adjusts for maximum mobility. The **Action³** manoeuvres smoothly due to its light rigid and well-balanced frame, and an improved ergonomic design that enables the user to become more independent.





➤ **Improved function**

Emphasis on user and therapist needs has resulted in easier and finer adjustments within the frame. These increased configuration choices ensure improved postural position maximising user function and independence.

➤ **Performance**

The *Invacare*® **Action³** is designed to ensure optimal driving performance with the lightweight rigid frame and high quality components providing easy rolling and manoeuvrability.

➤ **Functional**

Adaptable and modular, the **Action³** is ideally suited for accommodating progressive changes in user needs. It is also simple to refurbish with a range of options.

➤ **Easily transportable**

Simple to fold and dismantle with a low transfer weight, the **Action³** can be transported in a small car.





14,5 kg weight

The total weight of 14,5 kg addresses the need for improved handling, storage and transportation. The transportation weight (without legrest and rearwheel) is 9 Kg.



Brake and Brake extension

Brake application is improved by increased brake lever length. The brake handle can be easily folded for side transfers.



Length compensating Leg support

The new elevating legrest offers optimal length compensation with the elevation stepped in approx. 3 - 5 degrees. The weight allowance of the legrest is 60 kg in accordance to ISO standard.




NEW "quick release" Foot support

The swing in / swing out foot support is equipped with an easily operated, quick release mechanism. The height of the foot support can be adjusted in steps of 1 cm with two different inner tubes to meet low and high seat height requirements. Smooth new head support fixation with no flushing.



Active, standard, medium and amputee positions

Invacare® **Action³** offers as standard the possibility of increased driving performance with a choice of 20 different rear wheel positions. The rear wheel can be positioned in a standard or medium position, or in either the active or amputee position to suit the user need.

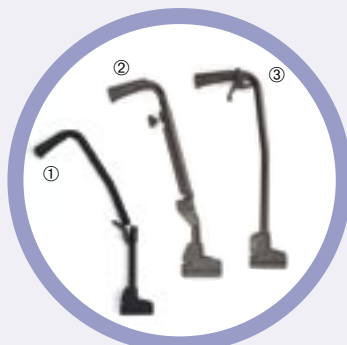

Backrest types

On top of the fixed backrest (0° or 7°), **Action³** offers 4 different types of backrests to suit user's needs:
 ① folding, ② angle adjustable (-6° / +24°) and ③ racket reclining (0°/30°).



New shape for a better support

The new head support provides ease of angle adjustment and a high support. Supplied as standard with a 3° seat angle, **Action³** seat angle can be adjusted to 6° (active user) or 0° (hemiplegic user) to facilitate foot propulsion.

Options



Hubbrake

Attendant and dual operation hubbrake available.



Transit chair

A transit version is available offering seat heights of 46, 48 and 51 cm.



Tension Adjustable Back support

Provides the user with increased comfort and appropriate spinal support.



Armrest options

Available in fixed height adjustable version, or swing-back version, variable from 18-28 cm with a short or long ampad.



Increased safety while seated

Buckled or velcro belts are available.



Cushions

Standard and gel cushions available.



Foot plates

The rigid foot supports are ISO approved. Available in 2 versions; fixed foot support or an angle and depth adjustable option, with 5-cm depth infinite angle adjustment.



Front castor

Front castors available in 6" and 8" pneumatic or solid tyre, various options of tyre width also available.



Rear wheel

Rear wheels are available in 22" or 24" composite, spoked, solid and pneumatic.

Options



Handrim

Choice of plastic coated, anodized or steel coated.



Amputee legrest

Height, angle and depth adjustable, the amputee legrest is also swing-away and should be associated with anti-tippers.



Anti-tip castors

Provide extra stability on slopes or uneven ground.



Transparent tray

Useful for reading or other activities, including mealtimes.



Plaster legrest

The *Invacare® Action³* can be equipped with plaster legrest and combined to anti-tippers. It offers then a fixed and stable 90° support.



Cane holder



Hemiplegic support


The *Invacare® Action³* can be fitted with dual handrim, low seat height, brake extension and hemiplegic arm support.



One arm lever drive

Control system allowing the user with low muscular tone to propel the wheelchair. Using one arm through the lever pendular motion which propels it forward or backwards, steering and braking integrated.

Technical data

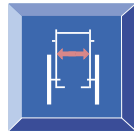


38, 40.5, 43, 45.5,
48, 50.5 cm




40, 45 cm




41, 43.5, 46, 48.5,
51 cm

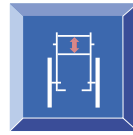

40, 43, 51 cm

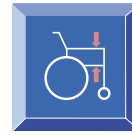

24 / 18 - 28 cm

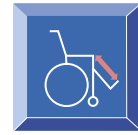


35 - 48 cm
34 - 47 cm




80° / 0° - 80°



0° / 7° / 0° - 30° /
-6° - +24°

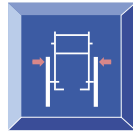

57 - 69.5 cm

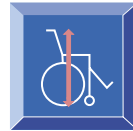


81.5 - 102 cm




100 cm



85 cm

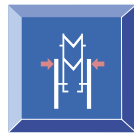

28.5 cm

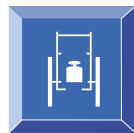

14.5 kg



max. 125 kg

Frame colour



Sapphire Blue



Anthracite Black



Dark Red



Pearl Grey

Please note colours may vary slightly from those displayed above.

Upholstery



Black nylon

All information quoted is believed to be correct at time of print.
Invacare® reserves the right to alter product specifications without prior consultation.



INVACARE® Limited
South Road · Bridgend Industrial Estate · BRIDGEND CF31 3PY
Tel. +44 (0) 1656 647327 · Fax +44 (0) 1656 649016
E-mail: uk@invacare.com · www.invacare.co.uk

