

QUANTUM[®] 6000



Q6000
Shown in Red with
Euro Comfort seat and
PG Remote Plus controller

QUANTUM[®] 6000

The Quantum[®] 6000 incorporates aggressive new features, powerful 4-pole Hammer motors, and our patent pending Mid-Wheel 6™ design platform to deliver unmatched performance and tight-space manoeuvrability while delivering limitless rehab capability.



Q6000
Shown in Red with
Euro Comfort seat and
PG Remote Plus controller

Pride
Mobility Products Ltd.
www.pridemobility.com

Unit 106, Heyford Park, Camp Road
Upper Heyford, Oxfordshire, OX25 5HA
Tel: 44 1869 233 100
Fax: 44 1869 233 400



Specifications:

WEIGHT CAPACITY:	136kg (21 STONE)
MAXIMUM SPEED:	Up to 10.5 km/h (6.5 mph)
GROUND CLEARANCE:	9 cm (3.5")
TURNING RADIUS:	56 cm (22") <i>Without Foot Riggings</i>
LENGTH:	93.5 cm (38") <i>Without Foot Riggings</i>
WIDTH:	65 cm (25.75")
DRIVE WHEELS:	35.5 cm (14") pneumatic
FRONT CASTORS:	13 cm (5") solid
REAR CASTORS:	15 cm (6") solid
SUSPENSION:	Active-Trac® ATX Suspension with integrated rear dual-strut suspension and articulating castor beam
DRIVE TRAIN:	2-motor, mid-wheel drive
BRAKING SYSTEM:	Intelligent braking (electronic, regenerative disc brakes)
ELECTRONICS:	100 amp, PG Remote Plus controller
STANDARD SEATING:	Modular Seat
STANDARD FOOTREST:	Swing-away legrests
BATTERY CHARGER	Off-board, 8 amp (standard)
PER-CHARGE RANGE:**	Up to 40km (25 miles)
BATTERIES:	Two 12 volt, 70 amp
BATTERY WEIGHT:	24kg (53.5lbs.) each
BASE WEIGHT:	71 kg (157.5lbs.)
SEAT WEIGHT:	17 kg (37lbs.) Modular Seat
MAX. SEATING	
DIMENSIONS:	See prescription form for seat sizing
WARRANTY:	Lifetime warranty on frame 1-year limited warranty on electronics 2-year limited warranty on drive train 18-month PG Drive control systems



■ Mid-Wheel 6 allows six wheels on the ground



■ OMNI-Casters on front and rear



■ New Hammer Motors

Features:



- Mid-Wheel 6 allows six wheels on the ground for maximum stability
- Compatible with TRU-Balance® Power Positioning
- OMNI-Casters (nylon spherical-shaped casters) on front and rear prevent wheel hang-ups
- Side-mounted, easily accessible freewheel levers
- ATX Suspension (Active-Trac with extra stability) incorporates front OMNI-Casters for enhanced performance over more varied terrain
- 100 amp, PG Remote Plus controller
- 8 amp, off-board charger
- 4-pole Hammer motors
- Modular Seat (standard)

Options:

- Solid tyres
- Swing-away joystick
- Multi-axis joystick mount
- Elevating leg rests
- Oxygen tank holder
- Neoprene leg rest sleeves
- Vent tray
- TRU-Balance® power positioning
- TRU-Balance® power elevating seat
- Angle-adjustable footplates
- Lap belt
- Contoured seating
- High-mount foot platform
- 70° HD drop-in footrests
- Residual limb supports
- Contoured seating with solid seat pan


■ Swing-away joystick



■ Elevating leg rests



Colours:

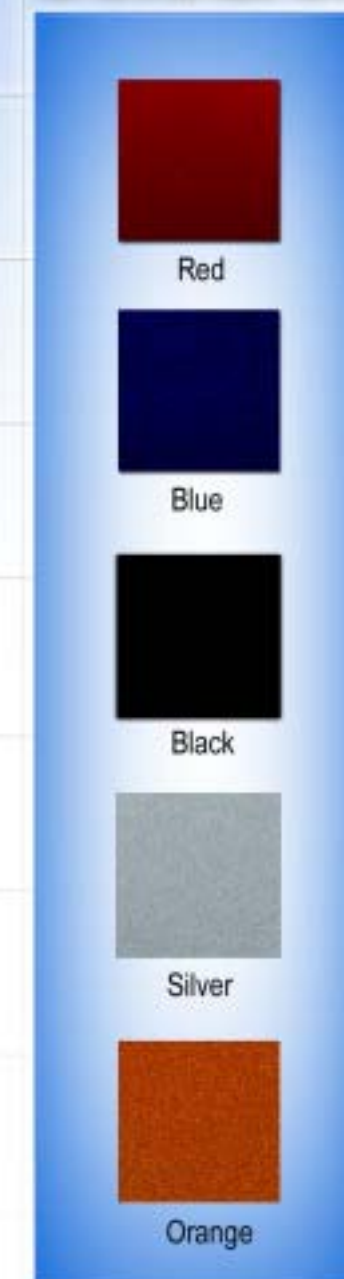

■ Oxygen tank holder



**The information contained herein is correct at the time of publication; due to the manufacturer's commitment to constant improvement and development, they reserve the right to alter specifications without prior notice.

***Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressures, ambient temperature, and the gradient and surface of the road or pavement.

***Exact product configurations shown on this literature may not be available in all markets